# Digital Graphic Printer

ご使用になる前に	JP
Before Using this Printer	GB
Avant d'utiliser cette imprimante	FR
Vor Verwendung dieses Druckers	DE
Prima di usare la stampante	IT
Antes de utilizar la impresora	ES
Alvorens de printer te gebruiken	NL
Innan du använder skri varen	SE
Antes de usar a impressora	PT
Ennen tulostimen käyttöä	FI
Før printeren tages i brug	DK
Πριν Χρησιμοποιήσετε τον Εκτυπωτή	GR
Før du bruker skr iveren	NO

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# WARNING

To prevent fire or shock hazard, do not expose the unit to rain or moisture.

To avoid electrical shock, do not open the cabinet. Refer servicing to qualified personnel only.

THIS APPARATUS MUST BE EARTHED.

**Symbol on the products**

This symbol indicates the equipotential terminal which brings the various parts of a system to the same potential.

This symbol is intended to alert the user to the presence of important operating and maintenance (servicing) instructions in the literature accompanying the appliance.

## For the customers in the U.S.A.

This equipment has been tested and found to comply with the limits for a Class A digital device, pursuant to Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference when the equipment is operated in a commercial environment. This equipment generates, uses, and can radiate radio frequency energy and, if not installed and used in accordance with the instruction manual, may cause harmful interference to radio communications. Operation of this equipment in a residential area is likely to cause harmful interference in which case the user will be required to correct the interference at his own expense.

You are cautioned that any changes or modifications not expressly approved in this manual could void your authority to operate this equipment.

This device requires shielded interface cables to comply with FCC emission limits.

## For the customers in Canada


This unit has been certified according to Standard CSA C22.2 NO.601.1.

## Important safeguards/notices for use in the medical environments

- 1 All the equipments connected to this unit shall be certified according to Standard IEC60601-1, IEC60950, IEC60065 or other IEC/ISO Standards applicable to the equipments.

- 2 When this unit is used together with other equipment in the patient area\*, the equipment shall be either powered by an isolation transformer or connected via an additional protective earth terminal to system ground unless it is certified according to Standard IEC60601-1.

\* Patient Area

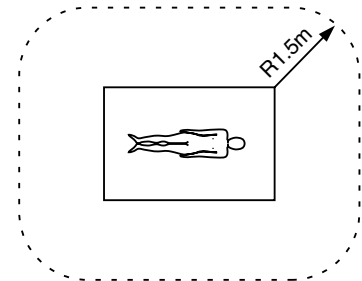- 3 The leakage current could increase when connected to other equipment.
- 4 This equipment generates, uses, and can radiate frequency energy. If it is not installed and used in accordance with the instruction manual, it may cause interference to other equipment. If this unit causes interference (which can be determined by unplugging the power cord from the unit), try these measures: Relocate the unit with respect to the susceptible equipment. Plug this unit and the susceptible equipment into different branch circuit. Consult your dealer. (According to Standard EN60601-1-2 and CISPR11, Class B, Group 1)

## Caution

When you dispose of the unit or accessories, you must obey the law in the relative area or country and the regulation in the relative hospital.

## Warning on power connection

Use a proper power cord for your local power supply.

## Warning on power connection for medical use

Please use the following power supply cord.

With connectors (plug or female) and cord types other than those indicated in this table, use the power supply cord that is approved for use in your area.

	United States	Canada
Plug type	HOSPITAL GRADE*	HOSPITAL GRADE*
Female end	E62405, E35708	LR53182, LL022442, LL088408
Cord type	E159216, E35496 Min.Type SJT Min.18 AWG	LL112007-1, LL20262, LL32121, LL84494 Min.Type SJT Min.18AWG
Minimum cord set rating	10A/125V	10A/125V
Safety approval	UL Listed	CSA

\*Note: Grounding reliability can only be achieved when the equipment is connected to an equivalent receptacle marked 'Hospital Only' or 'Hospital Grade'.

# Precautions

## On safety

- Check the operating voltage before operation. Operate the unit only with a power source specified in “Specifications”.
- Stop operation immediately if any liquid or solid object falls into the cabinet. Unplug the unit and have it checked by qualified personnel.
- Unplug the unit from the wall outlet if you will not be using it for a long time. Disconnect the power cord by grasping the plug. Never pull the cord itself.
- Do not disassemble the cabinet. Refer servicing to qualified personnel only.
- Do not touch the cutting blade of the printer.
- Connect the power plug of the printer to a wall outlet with a protective earth terminal. The safety earth should be properly established.

## On operation

Do not turn the power off while the printer is printing. The thermal head may be damaged.

## On the printer carriage

Do not carry or move the printer when the paper roll is placed in the printer. Doing so may cause a malfunction.


## On installation

- Place the printer on a level and stable surface. If you use the printer with placed on an uneven surface, malfunction of the printer is likely to occur.
- Do not install the printer near heat sources. Avoid locations near radiators or air ducts, or places subject to direct sunlight or excessive dust, humidity, mechanical shock or vibration.
- Provide adequate air circulation to prevent heat build-up. Do not place the printer on surfaces such as rugs, blankets, etc., or near materials such as curtains and draperies.
- If the printer is subjected to wide and sudden changes in temperature, such as when it is moved from a cold room to a warm room or when it is left in a room with a heater that tends to produce large amounts of moisture, condensation may form inside the printer. In such cases the printer will probably not work properly, and may even develop a fault if you persist in using it. If moisture or condensation forms, turn off the power and allow the printer to stand for at least one hour.

# Specifications

## Printer

Thermal head	Thin-film thermal head, 1280 dots
Gradation	256-level gradation (8-bit)
Memory capacity	4096 × 1280 × 8 (bit)
Picture size	320 × 100 mm (max.)
Printing speed	About 5.1 sec./image (1280 × 1280 dots)
Picture elements	4096 × 1280 dots (max.)
Dot density	H: 100 mm/1280 dots = 0.078 V: 100 mm/1280 dots = 0.078 Square dots

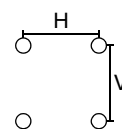Power supply	100 to 120V AC, 50/60Hz, 1.5A 220 to 240V AC, 50/60Hz, 0.8A
Temperatures	Operating 5°C ~ 35°C (41°F ~ 176°F) Storage/Transport -20°C ~ 60°C (-4°F ~ 140°F)
Humidity	Operating 20% ~ 80% Storage/Transport 20% ~ 90%
Dimensions	154 × 105 × 260 mm (w × h × d) 6.1 × 4.1 × 10.2 in. (w × h × d)
Weight	3.4 kg (1.54 lb.) (printer only)
Accessories provided	Printer paper UPP-110HG × 1 roll AC power cord × 1 Thermal head cleaning sheet × 1 Printer driver/Instructions for use disc × 1 Before Using this Printer × 1 Service Contact List × 1 Software License × 1 Media label × 1

GB

## Interface

### Parallel Interface

Control connector

IEEE1284-B Connector (D-sub 36-pin)

Input: 5V (TTL) (max.)

Output: 5V (TTL) (max.)

Data transmission

8-bit, parallel, IEEE STD 1284-1994 Standard

Compatible, Reverse nibble, ECP

Logic level

TTL

### Parallel Interface Connector Pin Assignment

Pin	I/O	Signal		
		Interface Mode		
		Compatible	Nibble	ECP
1	I	nStrobe	HostClk	HostClk
2	I/O	Data1 (LSB)		
3	I/O	Data2		
4	I/O	Data3		
5	I/O	Data4		
6	I/O	Data5		
7	I/O	Data6		
8	I/O	Data7		
9	I/O	Data8 (MSB)		
10	O	nACK	PtrClk	PeriphClk
11	O	Busy	PtrBusy	PeriphAck
12	O	Perror	AckDataReq	nAckReverse
13	O	Select	Xflag	Xflag
14	I	nAutoFd	HostBusy	HostAck
15		Not defined		
16-17		GND		
18	O	Peripheral Logic High (pull up to +5V at 1k $\Omega$ )		
19-30		GND		
31	I	nInit	nInit	nReverseRequest
32	O	nFault	nDataAvail	nPeriphRequest
33		Not defined		
34		Not defined		
35		Not defined		
36	I	nSelectIn	IEEE 1284 Active	IEEE 1284 Active

UP-D895 is compatible with the Compatible, Reverse Nibble, and ECP Modes with the bi-directional parallel interface (IEEE STD 1284-1994).

## USB Interface

Data transmission

Universal Serial Bus Specification Rev. 1.0 Standard

### USB Interface Connector Pin Assignment

Pin	I/O	Signal	Function
		VCC	Cable power, Max. current 100mA
	I/O	–Data	Data
	I/O	+Data	Data, pull up to +3.3V through a 1.5k $\Omega$ resistance
		Ground	Cable ground

The design and these specifications are subject to change without prior notification.

## Medical Specifications

Protection against electric shock:

Class I

Protection against harmful ingress of water:

Ordinary

Degree of safety in the presence of flammable anesthetics or oxygen:

Not suitable for use in the presence of flammable anesthetics or oxygen

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